


HYDRIM Wash Test Indicators

Monitoring Instrument Washers

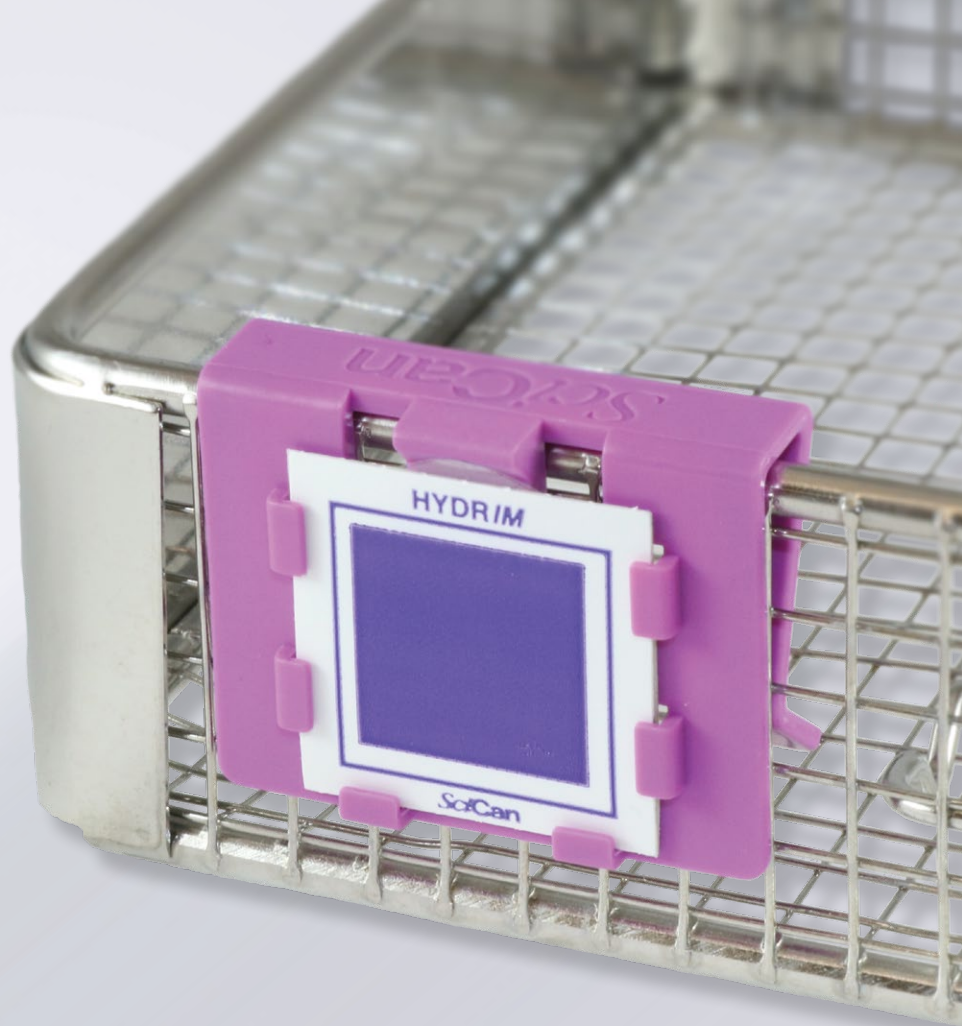

Think of it as a routine check up for your washer

The **HYDRIM Wash Test Indicators (WTIs)** have been designed to provide two-dimensional validation of HYDRIM washing efficacy. It is a simple, easy-to-use indicator that confirms instruments are effectively and consistently cleaned with HIP Cleaning Solution.



Safe

Synthetic test soil (non-toxic indicator) designed to monitor the efficacy of HYDRIM washers.



Reliable

Designed to provide verification of a consistently effective cleaning process in HYDRIM washers.



Smart Design

Indicators are folded “back-to-back” and placed inside the holder, monitoring spray performance from two directions.



Easy

WTI wash off characteristics clearly indicate whether the wash cycle has passed or failed; clean indicator at the end of the cycle indicates an effective cleaning process.

SciCan | Your Infection
Control Specialist

The only validated and effective cleaner in the market for use in HYDRIM units.

HIP Cleaning Solution is designed to remove any type of protein/debris from instruments. HIP is automatically added at various stages in the cycle to optimize the cleaning and protection of instruments. To maximize HYDRIM cleaning process, it is essential to use the certified and validated HIP Cleaning Solution that works for your HYDRIM unit during the wash cycle.



HIP Cleaning Solutions for HYDRIM 112W, L110W, C61W

HIP Solution is designed with the user in mind. Delivered in a format that fits inside the unit instead of taking up valuable cabinet space, the strong, built-in tube connector prevents leaks and allows for easy solution changes.

Part Numbers

PART NUMBER	DESCRIPTION	PART NUMBER	DESCRIPTION
WTI-SK	Wash Test Starter Kit (One Holder and 32 Indicators)	CS-HIPC-U	HIP Ultra Solution for C61W/WD G4 (8x750mL)
WTI-RK	Wash Test Refill Kit (96 Indicators)	G4-W10001	HYDRIM 112W G4+
CS-HIPT	HIP Ultra Solution for 112W G4+ (4x1.75L)	C61W-D01	HYDRIM C61W G4
CS-HIPL	HIP Cleaning Solution for L110W G4 (2x3.8L)	L110W-D01-G4	HYDRIM L110W G4

Coltene International Dental Group

Manufactured for:
Dent4You AG
Bahnhofstrasse 2
CH-9435 Heerbrugg

Distributed by:
SciCan Ltd.
1440 Don Mills Rd
Toronto, ON
M3B 3P9 Canada

Distributed by:
Coltène/Whaledent Inc.
235 Ascot Pkwy
Cuyahoga Falls, OH
44223, US

 **COLTENEGROUP**

HYDRIM is a registered trademark of SciCan Ltd.

SciCan | Your Infection Control Specialist