

STATCOOL

Direct-to-Drain for STATIM

Robust engineering that protects your drain, cycle after cycle.

STATCOOL is one of the smallest, most compact sterilizer exhaust cooling systems available, designed to fit easily in any under-counter space, vertically or horizontally.

It uses a simple, non-electric wax motor that activates instantly for immediate cooling. Its inline temperature monitoring continually checks water temperature making sure your plumbing is always protected.



Protects plumbing by reducing drain exposure to high heat.



Eliminates steam and moisture buildup caused by storing hot waste bottles in enclosed cabinets.



Improves safety for staff by removing the need to handle hot, exhaust waste bottles.



Installs neatly inside cabinetry, fitting securely along the cabinet wall beneath the counter.



Cools exhaust wastewater for up to two STATIM® autoclaves.



Saves time with its non electric, self regulating design that has no bottles to empty and needs virtually no maintenance.



Compact size fits under any counter

The innovative STATCOOL is small enough to fit in the palm of your hand. Its tiny stature gives you a wide range of installation options. Using the mounting clips, it can be configured for vertical or horizontal installation in any steri-center.



2-Year Warranty

STATCOOL's durability and long-term performance is backed by a 2-year warranty that protects your investment.



Digital Thermal Alarm

An inline water temperature monitor detects hot exhaust water ensuring STATCOOL always operates optimally.



Code-compliant installation

Installation hardware includes QAI-approved VersaCheck® QC dual check backflow preventer.

Water Systems Products & Accessories



STATPURE Water Filtration

Produces two types of high purity process water for sterilizers, automated instrument washers, dental bottles, and ultrasonic cleaners.

SWP400




DrainDock®

A multi-port drain adapter to create a unified drain connection point for steri-center equipment.

DD158 DrainDock 8-Port

DD158FK DrainDock 8-Port with fittings

DDFK1 DrainDock Fittings Kit

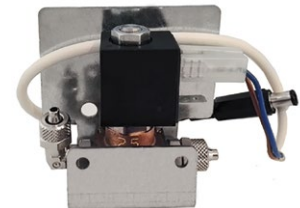

BRAVO G4 Auxiliary EV Kit

Replenish BRAVO G4's internal water reservoir with autoclave-quality water automatically — enabling truly hands-free convenience!

M7C200012

Distributed by:
SciCan Ltd.
1440 Don Mills Rd
Toronto, ON
M3B 3P9, Canada

Distributed by:
Coltene/Whaledent Inc.
235 Ascot Pkwy.
Cuyahoga Falls, OH
44223, USA

Manufactured for:
Dent4You AG
Bahnhofstrasse 2
CH-9435 Heerbrugg

For more information, visit www.scican.com

COLTENEGROUP

STATCOOL and STATPURE are trademarks of Coltene/Whaledent AG. DrainDock and VersaCheck are registered trademarks of Vista Water Group. STATIM is a registered trademark; BRAVO is a trademark of SciCan Ltd. This system incorporates proprietary technology and is protected by pending patent applications and associated intellectual property rights, and is provided under exclusive license.

SciCan | Your Infection
Control Specialist