

STATPURE

Water Filtration System

Simplified water filtration and easy maintenance to match your workflow.

STATPURE's unique four-filter configuration produces exceptional quality water for autoclaves while also making hyper-filtered RO water for dental bottles, ultrasonic cleaners, and final rinses in automated instrument washers.

The streamlined design neatly organizes tubing for easy, intuitive installation, while QuickTwist filters cut maintenance time down to moments.



ERGONOMIC

Autoclave Fill Wand

The wand's grip-pressure control reduces hand strain making daily autoclave filling easier and faster for your busy practice.



MORE TIME

STATPURE eliminates storing and handling heavy bottled water. Its simple design cuts filter change times in half, and it can be configured to autofill SciCan reprocessing equipment.



MORE SAVINGS

STATPURE lowers your clinic's costs long-term by eliminating bottled water purchases. Its four-filter system produces high-purity water without extra filters.



MORE WATER

Powered by your municipal water supply, STATPURE delivers a high volume of on demand water. Match your clinic's needs with the 4.25-gallon or 9.5-gallon net storage tanks to keep your instrument reprocessing moving.



2-Year Warranty

STATPURE is backed by a 2 year warranty that protects your investment and provides confidence you can rely on dependable water filtration for your Stericenter workflow.



Maintenance in a moment


STATPURE's QuickTwist filter design simplifies filter changes by eliminating the need to disconnect and reconnect tubing and fittings. Instead, simply unscrew the filter and replace.

Average 5-Year Filter Replacement*

Filter	Part No.	12 Mo	24 Mo	36 Mo	48 Mo	60 Mo
Prefilter	1 SWP01F	~	~	~	~	~
Hyperfiltration	2 SWP02F				~	
Deionization	3 SWP03F	~	~	~	~	~
Polishing	4 SWP04F	~	~	~	~	~

*Actual filter lifespan depends on local water quality, inbound water pressure, and volume of water used. Filter replacement intervals will be based on TDS readings.

Installation Space Requirements

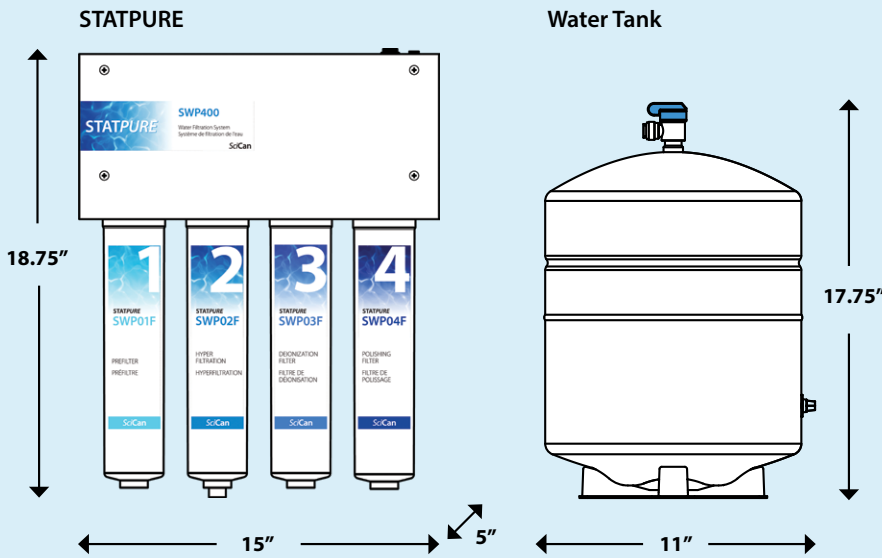

4.25 Gallon Net Water Tank STANDARD

11" Diameter
17.75" Height

SWP-TA425

9.5 Gallon Net Water Tank OPTIONAL

15.75" Diameter
22" Height

SWP-TA950

Filter Kits & Accessories

STATPURE Filter Kit 1
(Prefilter & Polishing)

SWP05KIT

STATPURE Filter Kit 2
(DI 2-Pack)

SWP06KIT

STATPURE Multipack Kit 3
(1 x Pre, 2 x DI, 1 x Pol)

SWP07KIT

Auxiliary EV Kit
for BRAVO G4

M7C200012

Distributed by:
SciCan Ltd.
1440 Don Mills Rd
Toronto, ON
M3B 3P9, Canada

Distributed by:
Coltene/Whaledent Inc.
235 Ascot Pkwy.
Cuyahoga Falls, OH
44223, USA

Manufactured for:
Dent4You AG
Bahnhofstrasse 2
CH-9435 Heerbrugg

For more information, visit www.scican.com

SciCan | Your Infection
Control Specialist

COLTENEGROUP

STATPURE is a trademark of Coltene/Whaledent AG. BRAVO is a trademark of SciCan Ltd. This system incorporates proprietary technology and is protected by pending patent applications and associated intellectual property rights, and is provided under exclusive license.