

VistaCool

Direct-to-Drain Thermal
Reduction System
for Autoclave Wastewater



Eliminates the need for autoclave condenser bottles

by cooling sterilizer wastewater and sending it directly down the drain automatically.



Protects Cabinets

Helps prevent mold, mildew, rot and rust caused by spills and humid environments created by exhaust from condenser bottles.



Protects Drains

Prevents autoclave wastewater from damaging plumbing. Built-in backflow prevention meets or exceeds local plumbing codes.



Protects Staff

Eliminates the need for office staff to empty hot water bottles manually.




Saves Time

Non-electric, self-regulating system saves staff time—there's no user maintenance required!



Install anywhere.

The innovative design of the new VistaCool 3 provides all the same benefits as its predecessors, while occupying less space. The system's included mounting clips and pressurized drain line enable it to be installed vertically or horizontally in almost any location.

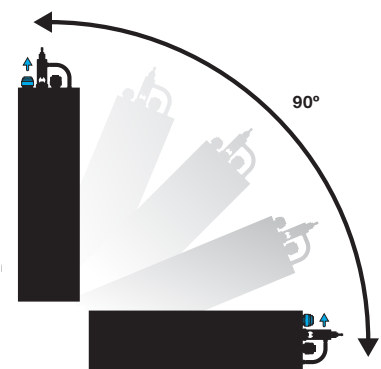

System Specifications



Note: Optional, adjustable tank stands raise the installed height of V7501 & V7502 models between 1/2" and 5" in half-inch increments. or S7545 for V7501; order S7547 for V7502.


Installation Requirements

1/4" connection to cold water supply: Supply fittings, water supply valve, cold water feed line and QAI-approved VistaCheck Dual Check Backflow Preventer(s) are all included! Air gap accessory available separately if required by local code authorities.

Drain connection with trap: Inline thermal sensor(s) included to ensure wastewater is cool enough to send directly down the drain. 1-1/2" utility tee with push-in fittings included with V7501 & V7502 models.

Optional DrainHub multi-port drain adapter recommended to allow multiple drain connections in one small, space-saving hub assembly.

Recommended Accessories

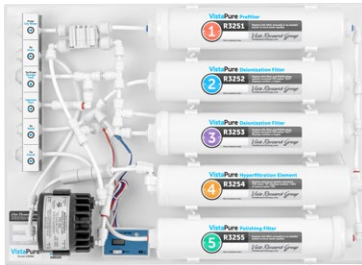

VistaPure

Water Filtration

Provides two grades of high-quality water for autoclaves, dental bottles, ultrasonic cleaners and the rinse cycle of instrument washers.

V3000



VistaPure AutoFill

Replenish any new or older-model STATIM autoclave with distilled-quality water automatically, enabling truly hands-free convenience.

S3020




DrainHub

Multi-Port Drain Adapter

This QAI-approved, patent-pending device accepts up to eight drain connections in one small hub—saving space, cost and time!

DH154 - Four-port drain adapter

DH158-82 - Eight port drain adapter with accessories kit

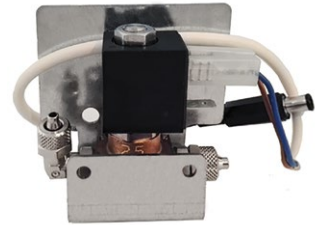

BRAVO G4

Auxiliary EV Kit

Replenish the SciCan BRAVO G4's internal water reservoir with autoclave-quality water automatically — enabling truly hands-free convenience!

M7C200012

Distributed by:
SciCan Ltd.
1440 Don Mills Rd
Toronto, ON
M3B 3P9, Canada
TF +1 (800) 667-7733
T +1 (416) 445-1600
customerservice@scican.com

Distributed by:
Coltene/Whaledent Inc.
235 Ascot Pkwy.
Cuyahoga Falls, OH
44223, USA

VistaPure and VistaCool are
manufactured for:
Crosstex International, Inc.
1625 Pinnacle Peak Rd
Phoenix, AZ 85027, USA
T +1 888-276-7783
HuFriedyGroup.com

For more information, visit www.scican.com

SciCan | Your Infection
Control Specialist

STATIM is a registered trademark and BRAVO is a trademark of SciCan Ltd. All Vista product names are trademarks of HuFriedy Mfg. Co. LLC, its affiliates or related companies, unless otherwise noted. Marks not registered in all jurisdictions