

SciCan



# VistaCool 3™

Direct-to-Drain Thermal Reduction System for Autoclave Wastewater

**Eliminates the need for autoclave condenser bottles by cooling sterilizer wastewater and sending it directly down the drain automatically.**



## Protects Cabinets

Helps prevent mold, mildew, rot and rust caused by spills and humid environments created by exhaust from condenser bottles.



## Protects Staff

Eliminates the need for office staff to empty hot water bottles manually.



## Protects Drains

Prevents autoclave wastewater from damaging plumbing. Built-in backflow prevention meets or exceeds local plumbing codes.




## Saves Time

Non-electric, self-regulating system saves staff time—there's no user maintenance required!




## Install anywhere.

The innovative design of the new VistaCool 3 provides all the same benefits as its predecessors, while occupying less space. The system's included mounting clips and pressurized drain line enable it to be installed vertically or horizontally in almost any location.

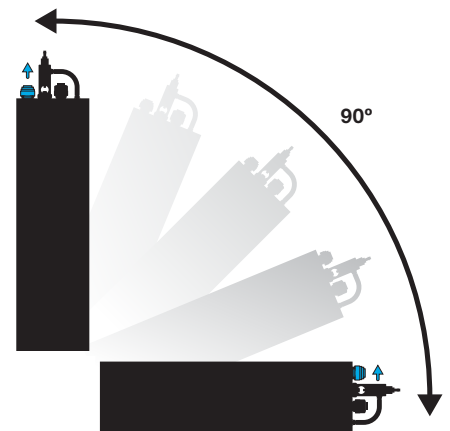

 **COLTENE**

## System Specifications

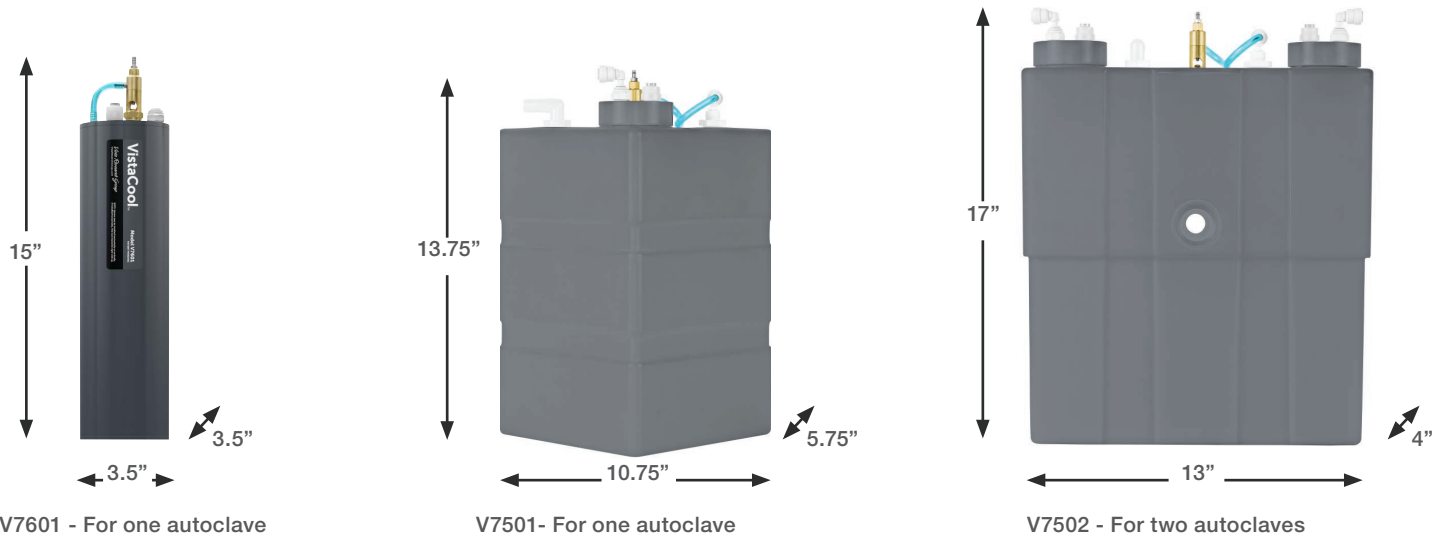

V7601 - For one autoclave

V7501- For one autoclave


V7502 - For two autoclaves

Note: Optional, adjustable tank stands raise the installed height of V7501 & V7502 models between 1/2" and 5" in half-inch increments. Order S7545 for V7501; order S7547 for V7502.

## Installation Requirements

**1/4" connection to cold water supply:** Supply fittings, water supply valve, cold water feed line and CSA-certified VistaCheck™ Dual Check Backflow Preventer(s) are all included! Air gap accessory available separately if required by local code authorities. Inline thermal sensor(s) included to ensure wastewater is cool enough to send directly down the drain. To be connected to WaterMark non-return valve according to local authorities (not supplied).

## Recommended Accessories

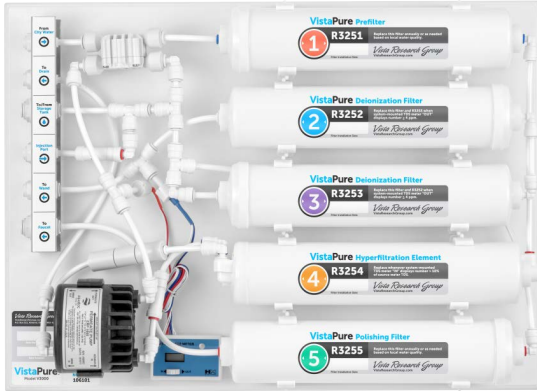

### VistaPure™ Water Filtration

Provides two grades of high-quality water for autoclaves, dental bottles, ultrasonic cleaners and the rinse cycle of instrument washers.

V3000

### VistaPure™ AutoFill

Replenish any new or older-model STATIM autoclave with distilled-quality water automatically, enabling truly hands-free convenience.

S3020 - AutoFill for STATIM

Distributed by:  
**SciCan Ltd.**  
1440 Don Mills Road  
North York, ON M3B 3P9 / Canada  
T +1 416 445 1600  
TF +1 800 667 7733  
customerservice@scican.com

VistaPure and VistaCool are manufactured by:  
**Crosstex International, Inc.**  
1733 W. Parkside Lane  
Phoenix, AZ 85027 / USA  
TF +1 800 531 2221  
VistaResearchGroup.com

[www.scican.com](http://www.scican.com)

 **COLTENE**

STATIM is a registered trademark of SciCan Ltd. VistaPure and VistaCool are trademarks of Vista Research Group.